

SENSICA

Sensilift Mini training

Sensilift Mini device

Wrinkle reduction






MY ANTI-AGING PREVENTION TOOL

- Suitable for all skin types
- Suitable for Men & Women
- Easy to use
- 3 energy levels
- Sensitive Touch ON/OFF button
- FDA cleared
- 2 temperature sensors
- Requires the use of Base Gel

Sensilift Mini

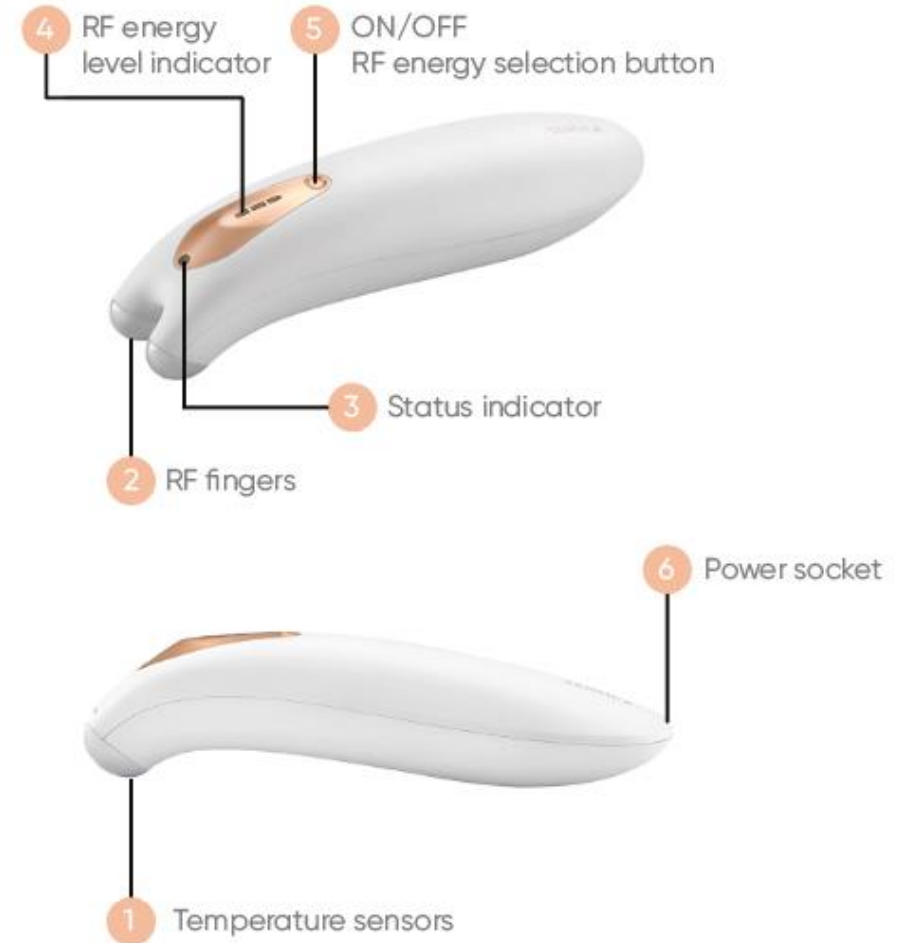

Benefits

- A professional skin tightening device
- Reduces the appearance of fine lines and wrinkles
- Skin feels tighter, firmer and renewed
- Stimulates collagen & elastin production
- Natural solution -Painless & non-invasive
- Featuring patented InTrense RFTM technology



How does RF (Radio Frequency) work?

- RF is a form of electromagnetic energy that heats the inner layer of skin while sparing the epidermis (outer skin layer)
- The RF waves penetrate deep into the skin to the source of collagen fibers where they produce a thermal effect. This thermal effect causes a remodeling of collagen fibers and in the long term produces new, and fresh collagen. This allows for long lasting results.
- New production of collagen is clinically proven!
- Additional effect is local improvements in blood circulation, thanks to the heat and the renewal of elastin fibers, further contributing to the skin's rejuvenation.

Safety

- Unique sensors constantly measure skin temperature and adjust delivery of RF energy accordingly
- RF delivery is deactivated if the treatment fingers are not in full contact with the skin



**Temperature and skin
contact sensors**

Base Gel

- The Base Gel is a light, cool and clear gel that enhances RF penetration through the skins layers, assisting treatment, and permitting effortless gliding over the skin
- The Base Gel is safe and easy to use
- Before starting treatment, clean and dry the skin area and apply an even layer of Base Gel
- After treatment, wipe off any residue with a damp cloth



⚠ The device cannot be used without the GEL!

Recommended treatment

Treatment areas:

Crow's feet (1)

Cheeks & around the mouth (2)

Chin & jawline (3)

Forehead (4)

Neck (5)

Décolletage & hands



- Sensilift Mini is intended for men and women who are generally healthy, and want to reduce wrinkles and tighten their skin
- Perform a preliminary test 48 hours prior to the first full treatment!
- Treat the desired areas with Sensilift Mini once a week for 8 consecutive weeks
- 5 minutes per area, 4 cm in size
- To sustain optimal long-term results, maintenance treatments may be performed as needed



Warning: Avoid treatment in the lower third of the neck, thyroid area, eyes cavities, nostrils and lips.

How to use

1.

Prepare your skin

Make sure the skin is clean, dry and free of any residue from lotions, creams or make-up



2.

Apply Base Gel

Apply an even layer of Base Gel



3.

Select energy level


Select the desired RF level



4.

Start treatment

Move the Sensilift Mini in slow circular motions over the treatment area

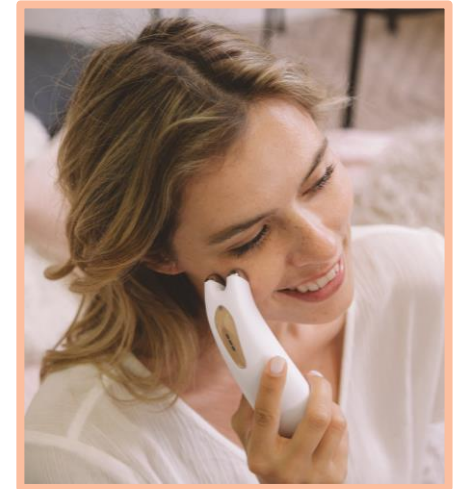

- Always begin with the lowest energy level and increase only if it feels comfortable
- Perform a preliminary test 48 hours prior to the first full treatment!
- Make sure that during the entire treatment the device remains in full contact with the skin while moving it in slow circular motions or simply back and forth



Thank you.